

CellVue®

Fluorescent Dye Kits for Cell Membrane Labeling

What are they?

CellVue® dyes are fluorescent probes for irreversible labeling of plasma membranes of live cells.


What do CellVue® Kits for Cell Membrane Labeling Offer?

- Versatility – use with any cell type or bioparticle with a membrane
- Provides stable labeling with minimal transfer from cell to cell
- Provides rapid, uniform membrane labeling
- Combine with fluorescent antibodies or markers of cell function
- Suitable for cell tracking and proliferation studies
- Several colors (UV to NIR) for multi-parameter studies (use with existing fluorochromes for more colors)
- Far-Red and NIR versions can provide greater signal to noise due to reduced background autofluorescence
- Compatible with flow cytometers, confocal and in vivo imaging equipment
- Convenient, easy-to-use kit format

How do they work?

The CellVue® cell linker kits use proprietary membrane labeling technology to stably incorporate a fluorescent dye with long aliphatic tails into the lipid regions of the cell membrane, see Figure 1. The labeling vehicle provided with the kit (Diluent C) is an iso-osmotic aqueous solution which contains no physiologic salts or buffers, detergents, or organic solvents and is designed to maintain cell viability while maximizing dye solubility and staining efficiency. The pattern of staining is dependent upon the cell type being labeled and the membranes of the cells.

Figure 1.

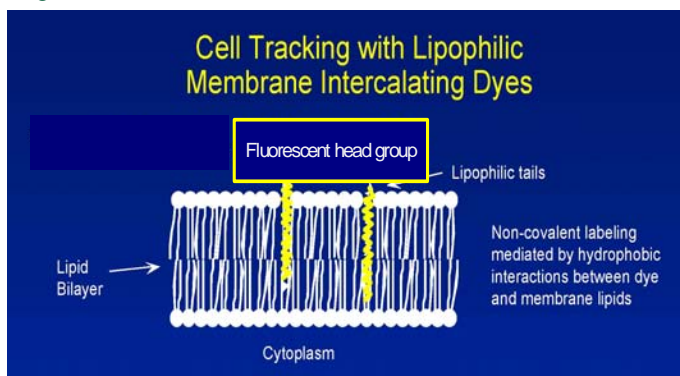


What are the advantages of Far Red and Near Infrared Fluorescence?

- Reduced autofluorescence background
- Greater signal to noise
- Greater ability to multiplex with visible probes due to minimal spectral overlap
- Excellent for use in combination with other cell tracking probes such as CFSE or PKH26

Methods for Cell Labeling with CellVue® Dyes

Mixing Steps and Critical Preparation Points

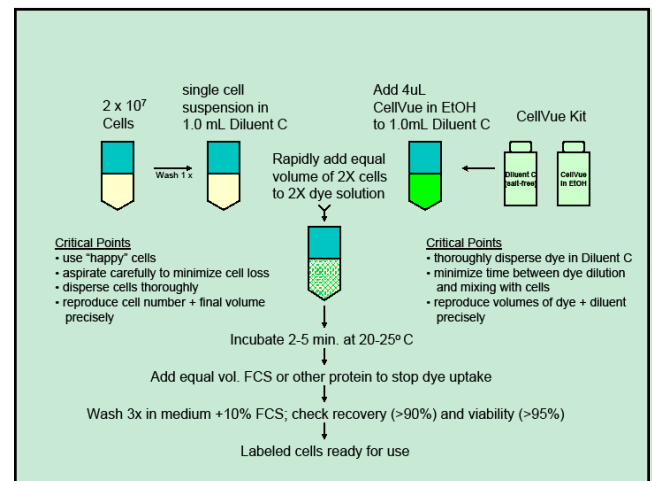

References:

- Gonzales et al. Mycobacterium tuberculosis H37 Rv induces ectosome release in human polymorphonuclear neutrophils. *Tuberculosis*, 90, 125-134 (2010).
- Katz SI et al. Sorafenib inhibits ERK1/2 and MCL-1 phosphorylation levels resulting in caspase-independent cell death in malignant pleural mesothelioma. *Can. Biol. & Ther.* 8:24, 101-111 (2009)
- Roy EJ et al. Imaging membrane intercalating near infrared dyes to track multiple cell populations. *J Immun. Meth.* 348, 18-29 (2009).
- Thomas DL et al. Experimental manipulations of afferent immune responses influence efferent immune responses to brain tumors. *Can. Immun. Immunother.* Sep. 57(9): 1323-33 (2008)
- Bantly AD et al. CellVue® Claret, a new Far-Red Dye, Facilitates Polychromatic Assessment of Immune Cell Proliferation. *Immun. Invest.*, 36, No 5-6, 581-605 (2007)
- Gertner-Dardenne J et al. Lipophilic Fluorochrome Trackers of Membrane Transfers between Immune Cells. *Immun. Invest.* 36, No 5-6, 665-685 (2007)
- Tario JD et al. Novel Lipophilic Tracking Dyes for Monitoring Cell Proliferation. *Immun. Invest.* 36, No 5-6, 861-885 (2007)
- Gertner J et al. Human gamma delta T lymphocytes strip and kill tumor cells simultaneously. *Immunology Letters*, 110(1): 42-53 (2007)
- Al-Mehdi AB et al. Increased Depth of Cellular Imaging in the Intact Lung using Far Red and Near Infrared Fluorescent Probes. *Int. J. Biomed Imaging*, 1-7 (2006).
- Stewart CC et al. Flow cytometer in the infrared: inexpensive modifications to a commercial instrument. *Cytometry*, Part A, 67A, #2, 104-111, (2005)

CellVue® Kits are available in 3 sizes:

Mini Kit (small): 0.1 ml of 1mM ethanolic dye stock and 1 x 10mL of diluent.

Midi Kit (medium): 0.2 ml of 1mM ethanolic dye stock and 6 x 10 ml of diluent.

Maxi Kit (large): 0.5 ml of 1mM ethanolic dye stock and 6 x 10 ml of diluent.

CELLVUE® PRODUCT LIST

Catalog Number	Name	Description	Price
C-1001	CellVue® Maroon	The CellVue Maroon Fluorescent Cell Linker Kit contains a 1 mM dye stock solution and cell labeling diluent. Dye fluorescence properties: Ex max = 647 nm and Em max = 667 nm. Provides stable labeling of the lipid regions of cell membranes.	Mini \$171.00 Midi \$359.00 Maxi \$563.00
C-1002	CellVue® Claret	The CellVue Claret Fluorescent Cell Linker Kit contains a 1 mM dye stock solution and cell labeling diluent. Dye fluorescence properties: Ex max = 655 nm and Em max = 675 nm. Provides stable labeling of the lipid regions of cell membranes.	Mini \$171.00 Midi \$359.00
C-1003	CellVue® Plum	The CellVue Plum Fluorescent Cell Linker Kit contains a 1 mM dye stock solution and cell labeling diluent. Dye fluorescence properties: Ex max = 652 nm and Em max = 671 nm. Provides stable labeling of the lipid regions of cell membranes.	Mini \$171.00 Midi \$359.00 Maxi \$563.00
C-1004	CellVue® Burgundy	The CellVue Burgundy Fluorescent Cell Linker Kit contains a 1 mM dye stock solution and cell labeling diluent. Dye fluorescence properties: Ex max = 683 nm and Em max = 707 nm. Provides stable labeling of the lipid regions of cell membranes.	Mini \$171.00 Midi \$359.00 Maxi \$563.00
C-1005	CellVue® Lavendar	The CellVue Lavender Fluorescent Cell Linker Kit contains a 1 mM dye stock solution and cell labeling diluent. Dye fluorescence properties: Ex max = 425 nm and Em max = 461 nm. Provides stable labeling of the lipid regions of cell membranes.	Mini \$171.00 Midi \$359.00 Maxi \$563.00
C-1006	CellVue® NIR815	The CellVue NIR815 Fluorescent Cell Linker Kit contains a 1 mM dye stock solution and cell labeling diluent. Dye fluorescence properties: Ex max = 786 nm and Em max = 814 nm. Provides stable labeling of the lipid regions of cell membranes.	Mini \$171.00 Midi \$359.00 Maxi \$563.00
C-1007	CellVue® NIR780	The CellVue NIR780 Fluorescent Cell Linker Kit contains a 1 mM dye stock solution and cell labeling diluent. Dye fluorescence properties: Abs max = 745 nm and Em max = 776 nm. Provides stable labeling of the lipid regions of cell membranes.	Mini \$171.00 Midi \$359.00 Maxi \$563.00
C-1008	Diluent C	Six vials containing 10 mL of Diluent C for membrane labeling with CellVue dyes.	\$154.00
C-1009	CellVue® Jade	The CellVue Jade Fluorescent Cell Linker Kit contains a 1 mM dye stock solution and cell labeling diluent. Dye fluorescence properties: Abs max = 478 nm and Em max = 508 nm. Provides stable labeling of the lipid regions of cell membranes.	Mini \$171.00 Midi \$359.00 Maxi \$563.00
C-1011	CellVue® Red	The CellVue Red Fluorescent Cell Linker Kit contains a 1 mM dye stock solution and cell labeling diluent. Dye fluorescence properties: Ex max = 567 nm and Em max = 588 nm. Provides stable labeling of the lipid regions of cell membranes.	Mini \$171.00 Midi \$359.00

CellVue® is a trademark of PTI Research Inc. used under license. CellVue® products are sold under sublicense from PTI Research, Inc.
US Patent Numbers 5,665,328; 7,462,347 B2 and 8,029,767 B2

COMPANY INFORMATION: Molecular Targeting Technologies, Inc., 833 Lincoln Ave., Unit 9, West Chester, PA 19380

P: 610.738.7938 F: 610.738.7928 website: www.mtarget.com



MOLECULAR TARGETING TECHNOLOGIES, INC.